





1. Open the first peel pouch and deliver the innermost sterile sealed pouch containing the allograft material to a sterile field.



2. Open the sterile sealed pouch and carefully deliver the allograft to the sterile field for transplantation.

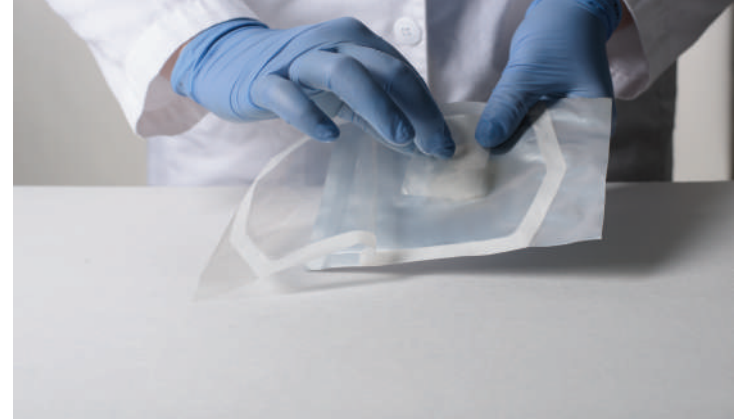

The allograft may be cut to a smaller size or shape with sterile scissors or scalpel if desired.

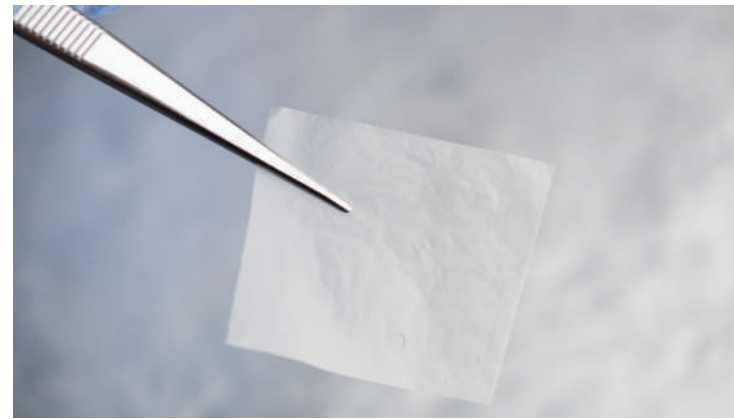

4. The allograft may be placed directly onto the surgical site and hydrated with sterile saline or other sterile solution, if desired, or hydrated for 1-2 seconds prior to application.

AmnioCore

AMNIOTIC MEMBRANE ALLOGRAFT

Product Offering

Type	Item Number	Size
Sheet	AMW-5230	2x3cm
	AMW-5330	3x3cm
	AMW-5340	3x4cm
	AMW-5440	4x4cm
	AMW-5460	4x6cm
	AMW-5480	4x8cm
	AMW-5660	6x6cm
	AMW-5690	6x9cm
	AMW-5616	6x16cm
	AMW-5920	9x20cm

The Stability Biologics® AmnioCore Product Line is regulated by the FDA under 21 CFR Part 1271 Human Cells, Tissues and Cellular and Tissue-Based Products (HCT/Ps). AmnioCore is processed by and donor eligibility determined by Stability Biologics®.

Stability Biologics® is registered with the FDA for tissue processing and accredited by the American Association of Tissue Banks (AATB).

1077 Central Parkway South Suite 500 San Antonio, TX 78232

© 2021 Stability Biologics® LLC. All Rights Reserved. www.stabilitybio.com
MKT-044 Rev. D Effective 09022021



Ordering Information

Phone: 855-267-5551

Fax: 855-267-4166

Email: cservice@stabilitybio.com